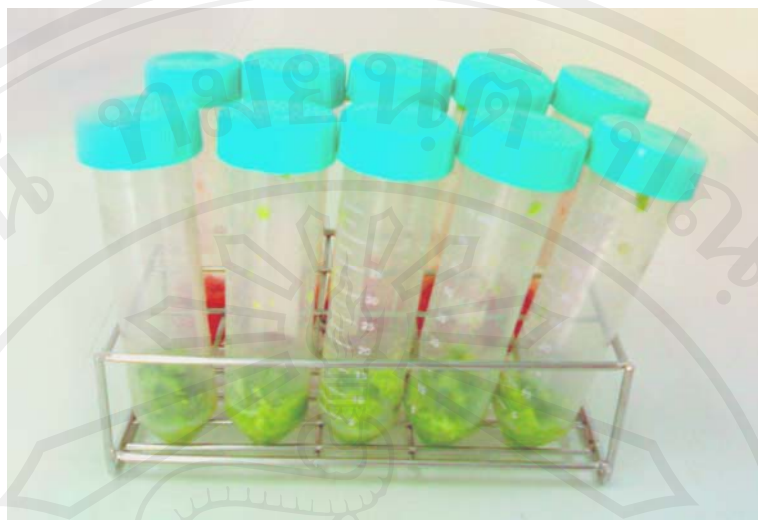




## APPENDIX

ลิขสิทธิ์มหาวิทยาลัยเชียงใหม่

Copyright© by Chiang Mai University  
All rights reserved



**Figure A-1** Weigh the sweet pepper samples at 10 gram/each and put them in 50-ml centrifuge tubes



**Figure A-2** Add 10 ml of acetonitrile and put them in 50-ml centrifuge tubes



**Figure A-3** Add 2 grams of sodium chloride for increase the efficiency of extraction



**Figure A-4** The extraction by homogenizer for extraction

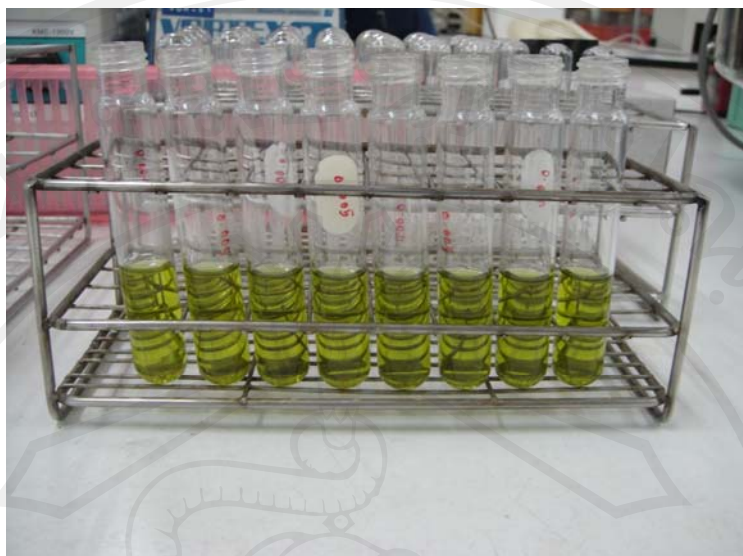


**Figure A-5** Centrifuge samples at the speed of 5000 rpm for 5 minutes



**Figure A-6** The cleaning up the samples by SPE alumina





**Figure A-7** The extracts for the clean up



**Figure A-8** Evaporate samples by nitrogen evaporator



**Figure A-9** The analyte samples by GC-MS

**CURRICULUM VITAE****NAME**

MRS.Juthamas Utrakian

**DATE OF BIRTH**

January 17, 1978

**EDUCATIONAL BACKGROUND**

B.Sc., Faculty of Sciences,  
Ratchaphat Chiang Mai University,  
Chiang Mai, Thailand (1996-2000)

**WORKING EXPERIENCE:**

2000 – Present

Scientist  
The Royal Project Foundation  
Chiang Mai, Thailand